

OXFORD CAMBRIDGE AND RSA EXAMINATIONS
A2 GCE
G495/01
PHYSICS B (ADVANCING PHYSICS)
Field and Particle Pictures
ADVANCE NOTICE

WEDNESDAY 23 JANUARY 2013: Afternoon
DURATION: 2 hours
plus your additional time allowance
MODIFIED ENLARGED 24pt

READ INSTRUCTIONS OVERLEAF

INSTRUCTIONS TO CANDIDATES

- **Take the article away and read through it carefully. Spend some time looking up any technical terms or phrases you do not understand. You are not required to research further the particular topic described in the article.**
- **For the examination on Wednesday 23 January 2013 you will be given a fresh copy of this article, together with a question paper. You will not be able to take your original copy into the examination with you.**
- **The values of standard physical constants will be given in the Advancing Physics Data, Formulae and Relationships booklet. Any additional data required are given in the appropriate question.**

INFORMATION FOR CANDIDATES

- **Questions in Section C of paper G495/01 Field and Particle Pictures will refer to this Advance Notice material and may give additional data related to it.**
- **Section C will be worth about 40 marks.**
- **Sections A and B of paper G495/01 will be worth about 60 marks.**
- **There will be 2 marks for quality of written communication (QWC) assessed in Sections B and C.**


THE LASER – A SURGEON’S BEST FRIEND

Since its invention and development over fifty years ago, the laser has found numerous applications in a wide range of contexts, including medicine. Although many different materials are used to make them, all lasers operate in essentially the same way. Electrons in the atoms of the selected material – called the gain medium – are excited to a higher energy level in such a way that the number of atoms with this energy E_2 is more than the number with a lower energy E_1 . This is called a population inversion: normally the number of atoms in the excited state is much smaller than the number in the lower energy state. If a photon of energy $hf = E_2 - E_1$ passes through the gain medium an electron in the higher energy level can be stimulated to fall down to the lower level. Another photon is emitted when this happens – see Fig. 1. The emitted photon will have exactly the same direction, frequency and phase as the photon which

5
10
15
20
25


stimulated its production in the first place: the two photons will be coherent. This process of stimulated emission was predicted by Einstein in 1917.

FIG. 1: STIMULATED EMISSION

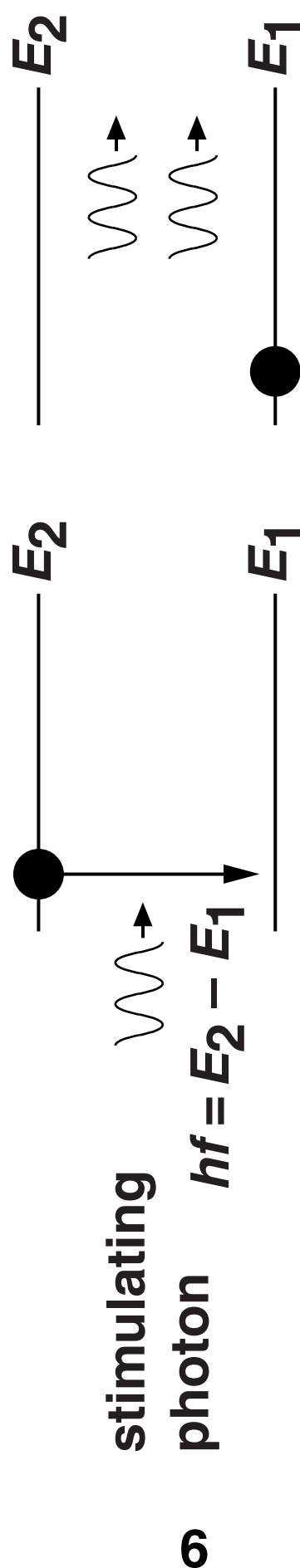

To make a laser, the gain medium is placed inside a 'cavity' made with two mirrors facing each other, Fig. 2. Photons travelling in a direction parallel to the cavity axis bounce back and forth within the cavity causing stimulated emission, increasing the number of photons all the time. A standing wave forms within the cavity which has a fixed length, d . One of the mirrors is a partial reflector so that light from the standing wave can leak through it and form the laser output. What emerges is a monochromatic, coherent and powerful beam of radiation.

FIG. 2: THE RADIATION REFLECTS BACK AND FORTH IN THE GAIN MEDIUM CAVITY, BUILDING UP A STANDING WAVE

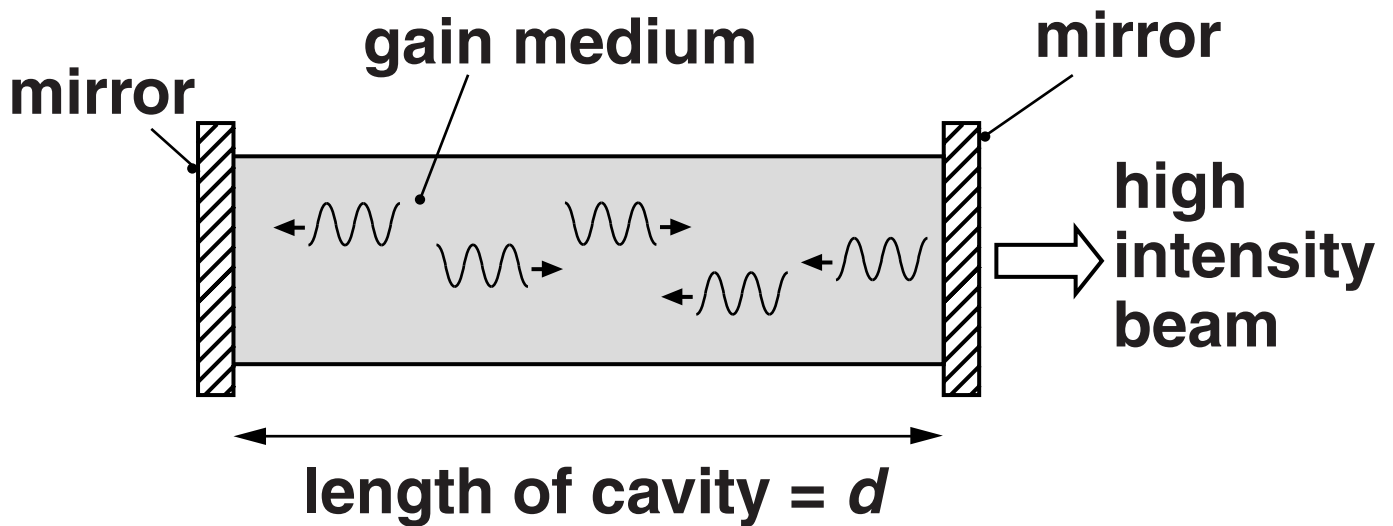

LASER ENERGY

45

In some lasers, the population inversion is maintained by providing the energy needed to raise electrons to the higher levels in a series of optical or electrical pulses. Thus the photons also emerge in regular bursts and this is called a pulsed laser. In others, the population inversion is maintained continuously and the laser radiation emerges as a steady beam – these are continuous wave (cw) lasers. The amount of energy deposited into the target from a laser depends upon the energy of the photons, the rate at which the photons arrive and the amount of time over which they are arriving.

50

55

60


LASERS IN SURGERY

The use to which a laser is put depends mostly on the wavelength of the radiation it emits. Human tissue responds in different ways to different wavelengths of the electromagnetic spectrum. Soft tissue is made of many different substances and each absorbs

65

light in its own way. Three substances 70
 in human skin and the soft tissue
 immediately beneath it are melanin
 (responsible for skin colouring), blood
 and water. The absorption by water
 of different wavelengths is shown in 75
 Fig. 3, in which the absorption
 coefficient, μ , is defined as the
 reciprocal of the depth within which the
 intensity falls by a factor of $1/e$ (about
 0.37). For example, an absorption 80
 coefficient of 0.1 cm^{-1} indicates that
 the intensity I of the radiation will
 have been reduced to about $0.37 I$
 within 10 cm. The absorption can be
 represented by the equation $I = I_0 e^{-\mu x}$, 85
 where I_0 is the incident laser intensity
 and x is the depth in cm.

FIG. 3: DIFFERENT WAVELENGTHS OF LASER RADIATION ARE ABSORBED BY DIFFERENT AMOUNTS BY WATER, ONE OF THE MAIN COMPONENTS (70%) OF SOFT TISSUE

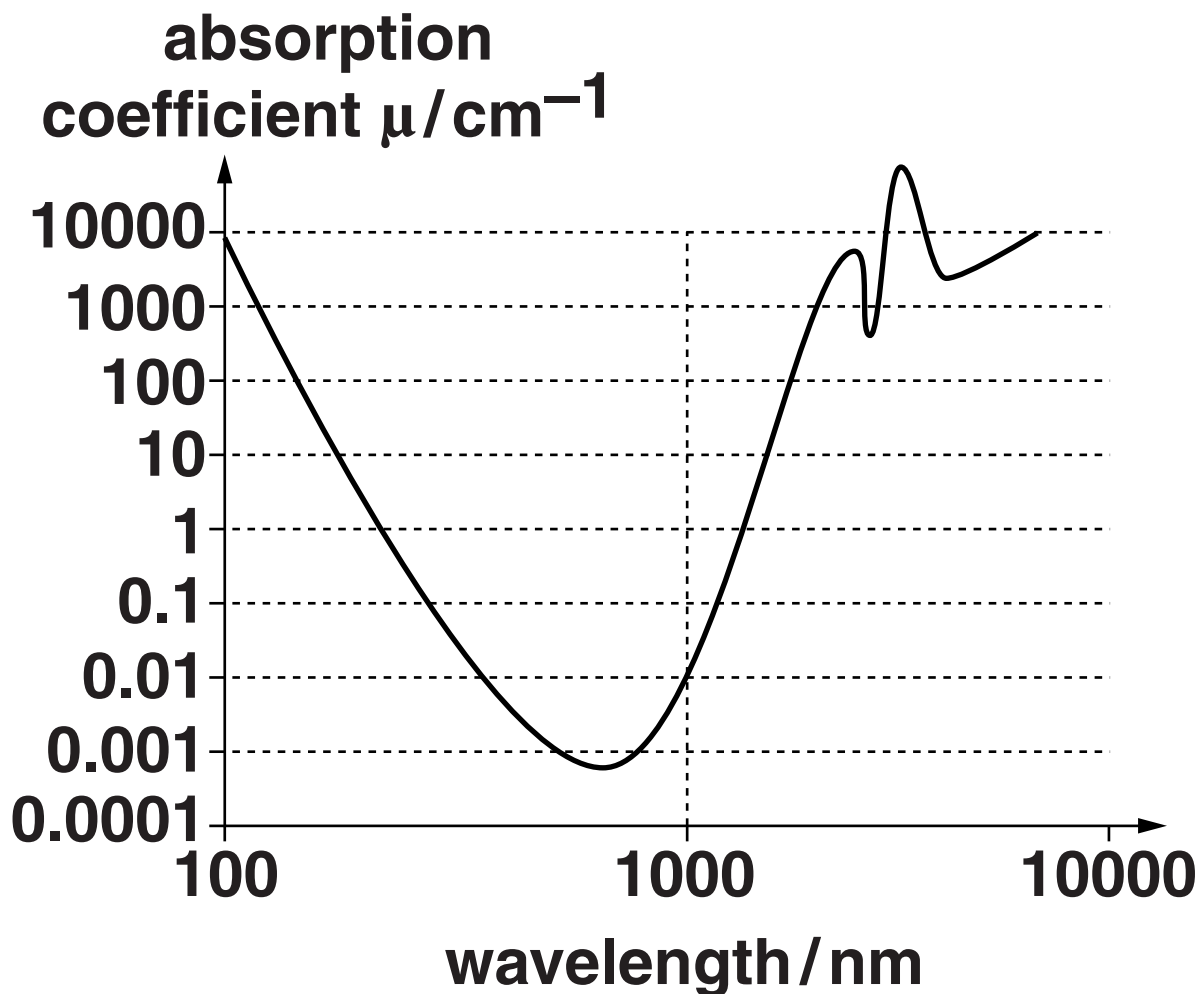

The interaction of laser radiation with such biological tissue occurs under the broad headings of thermal, mechanical, and chemical and the precise type of laser used in a surgical procedure will depend upon the particular type of surgery required. 90

THERMAL EFFECTS 95

For surgical procedures, tissue cells need to be broken up in a controlled fashion but without excessive heating taking place as this could cause damage to surrounding tissue. Heating tissue to about 60 °C causes what is known as coagulation, a process which causes it to shrink. This can be very useful, for example, for sealing blood vessels and preventing bleeding. A very narrow beam of powerful laser light can thus be used as a “bloodless scalpel”. At higher temperatures (greater than 100 °C), the water in tissue cells boils causing destruction of cells as the vapourised water expands, a process known as ablation. Carbon dioxide (CO₂) lasers, which 100 105 110

can be used in pulsed or continuous wave mode, have a wide range of applications, including straightforward surgical techniques involving coagulation. The CO₂ laser radiation is invisible, so the visible red light from a helium-neon (He-Ne) laser is used to provide a guide beam.

115

120

MECHANICAL EFFECTS

When particularly powerful laser beams are used, energy can be given to the tissue in a very short space of time. A pulse of high energy photons delivered in a nanosecond or so can break apart molecules without much heating taking place. This is called optical breakdown and leads to the production of short-lived “micro-plasmas”, tiny regions of ionisation with a density of liberated electrons as high as 10^{21} per cubic centimetre. The rapid expansion of the plasma leads to mechanical damage to the surrounding medium, although this remains a tiny region. Because it is a non-thermal effect and cell destruction


125

130

135

is highly localised this can be used 140
to remove tissue on a very delicate
scale, for example in the eye. One
type of laser used for this work is the
pulsed neodymium:YAG (Nd:YAG)
laser in which the gain medium is 145
an yttrium-aluminium-garnet crystal
doped with the element neodymium. Its
radiation is very penetrating and, when
administered in short pulses, is used
for eye surgery. 150

**FIG. 4: NEODYMIUM:YAG LASERS CAN
BE USED TO MAKE REPAIRS TO TINY
LOCALISED AREAS IN THE RETINA OF
THE EYE**

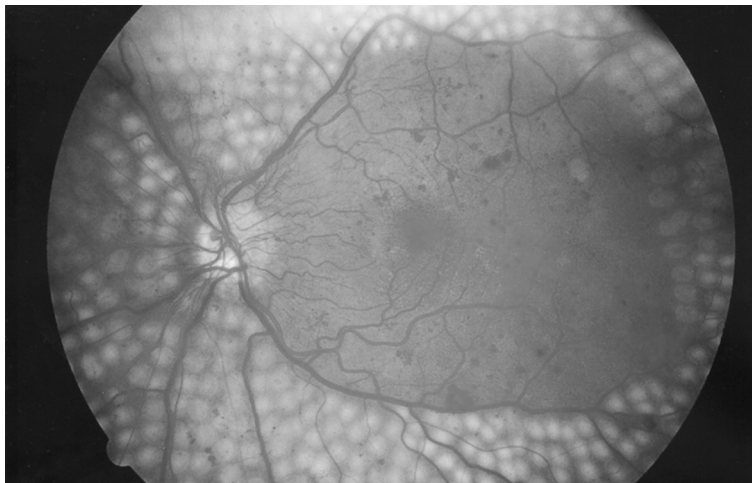

PHOTOCHEMICAL REACTIONS

The behaviour of some chemicals can be affected very profoundly by light and other parts of the electromagnetic spectrum. Certain drugs, for example, become very active when illuminated with short-wavelength visible light. If the chemical is reacting in something not far below the skin, then the laser can be shone directly onto the skin surface. However, to illuminate deeper regions, the light has to be directed using optical fibres. For this work, metal-vapour lasers are used, in which the gain medium consists of a mixture of neon and a vapourised metal, often gold.

155

160

165

THE FUTURE

By applying the same principles to different materials for more than fifty years, lasers have found uses in a very wide range of applications with many, such as DVD players, bar code readers and telecommunications, encountered every day. Current developments of

170

175

solid state lasers have started a new era in laser technology, though the basic principle remains the same. As higher-powered semiconductor lasers become a reality, no doubt the world of medicine and surgery will continue to benefit from this truly cutting-edge technology.

180

END OF ARTICLE



Copyright Information

OCR is committed to seeking permission to reproduce all third-party content that it uses in its assessment materials. OCR has attempted to identify and contact all copyright holders whose work is used in this paper. To avoid the issue of disclosure of answer-related information to candidates, all copyright acknowledgements are reproduced in the OCR Copyright Acknowledgements Booklet. This is produced for each series of examinations and is freely available to download from our public website (www.ocr.org.uk) after the live examination series.

If OCR has unwittingly failed to correctly acknowledge or clear any third-party content in this assessment material, OCR will be happy to correct its mistake at the earliest possible opportunity.

For queries or further information please contact the Copyright Team, First Floor, 9 Hills Road, Cambridge CB2 1GE.

OCR is part of the Cambridge Assessment Group; Cambridge Assessment is the brand name of University of Cambridge Local Examinations Syndicate (UCLES), which is itself a department of the University of Cambridge.

